

SIEMENS

MOBILETT Plus/ -E/ -M

SP

Planning Guide

© Siemens AG 1994

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

English

Doc. Gen. Date: 05.01

Print No.: RXR8-215.021.01.04.02

Replaces: RXR8-215.021.01.03.02

Chapter	Page	Revision
all	all	01
all	all	02
all	all	03
all	all	04

Document revision level

The document corresponds to the version/revision level effective at the time of system delivery.

Revisions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order current revision levels.

Disclaimer

The installation and service of equipment described herein is to be performed by qualified personnel who are employed by Siemens or one of its affiliates or who are otherwise authorized by Siemens or one of its affiliates to provide such services.


Assemblers and other persons who are not employed by or otherwise directly affiliated with or authorized by Siemens or one of its affiliates are directed to contact one of the local offices of Siemens or one of its affiliates before attempting installation or service procedures.


	Page
1 General Notes	1 - 1
General notes	1 - 1
Safety	1 - 2
System configuration	1 - 3
2 Room Planning	2 - 1
Dimensions: MOBILETT Plus/ -E/ -M	2 - 1
3 Technical Data	3 - 1
Electrical data	3 - 1
Weight and heat dissipation	3 - 1
Environmental conditions	3 - 1
Packing and transport routes	3 - 2
Surface colors	3 - 2
4 Changes to Previous Version	4 - 1

This page intentionally left blank.

General notes

- With distribution of these revision level, all preceding planning guides, Speed - Infos (PG's) and drafts lose their validity.
- All layouts issued by the Planning Departments must bear a note referring to the installation and delivery conditions of Siemens Medical Engineering Group. The installation and delivery conditions must be submitted with the layouts.
- Unless otherwise specified, all dimensions are indicated in "mm".

-  - The symbol indicates a change (see revision status).

-  - Orientation points
Points specific to system components to which reference is made when positioning system components to each other or in the room.
The isocenter of a radiographic system is always illustrated as the orientation point.

- Fixpoints
Clearly marked points on system components, installation ceiling, walls or floor on which cable outlets are located.
Illustration in the drawings: octagon with letter/number-combination.
The cable lengths specify the maximum fixpoint distances and thus the maximum distances between the individual system components.

- Room height
The room height is the distance measured from the top surface of the floor to the bottom surface of the ceiling structural elements (Unistrut rails) (bottom surface of drop ceiling).

- Room lighting
According to DIN 68 68-57 (international standard in preparation), the lighting in rooms in which image playback devices (monitors) are used for diagnosis, the following requirements must be met:

adjustable, no anti-glare screen, reproducible adjustment of the lighting (e. g. dimmer with scala),

no glare or reflection from windows, lights and light boxes in the standard working position of the monitors.


Hotline + 49 (9191) 18 - 8080

Safety ♦

- The provisions of the relevant fire safety regulations must be observed for the premises.
- The system has been developed according to EN 60601 - 1.
- Minimum dimensions (e. g. room heights, safety distances) indicated in the planning guides are marked "min."
- Basic strength against electromagnetic sources of interference.
Result of lightning discharges.
The protection targets of the different lightning protection areas up to the unit connection are also specified in the IEC 1024, DIN 48810, VDE 0675 and in the DEMVT recommendations.


System configuration

(schematic)

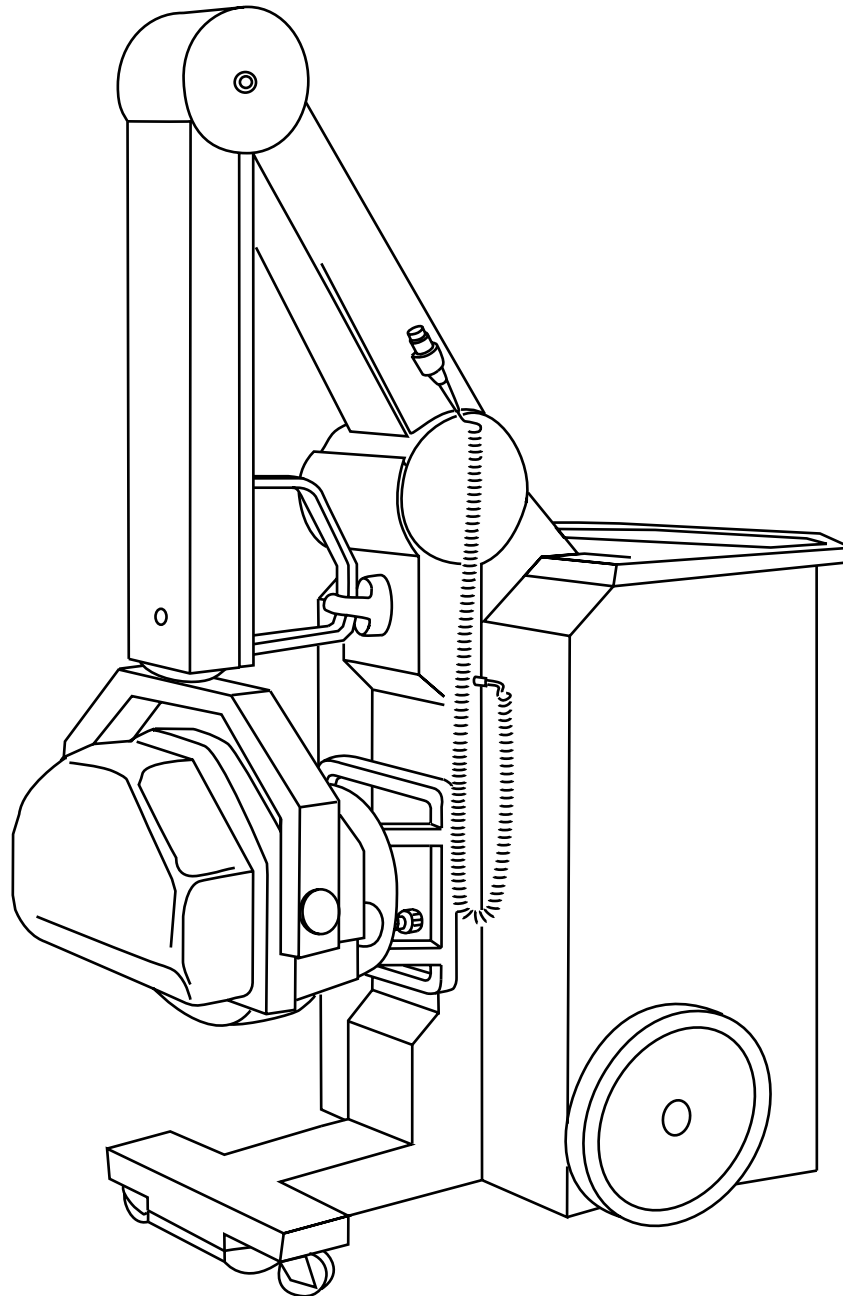

This page intentionally left blank.

Dimensions: MOBILETT Plus/ -E/ -M

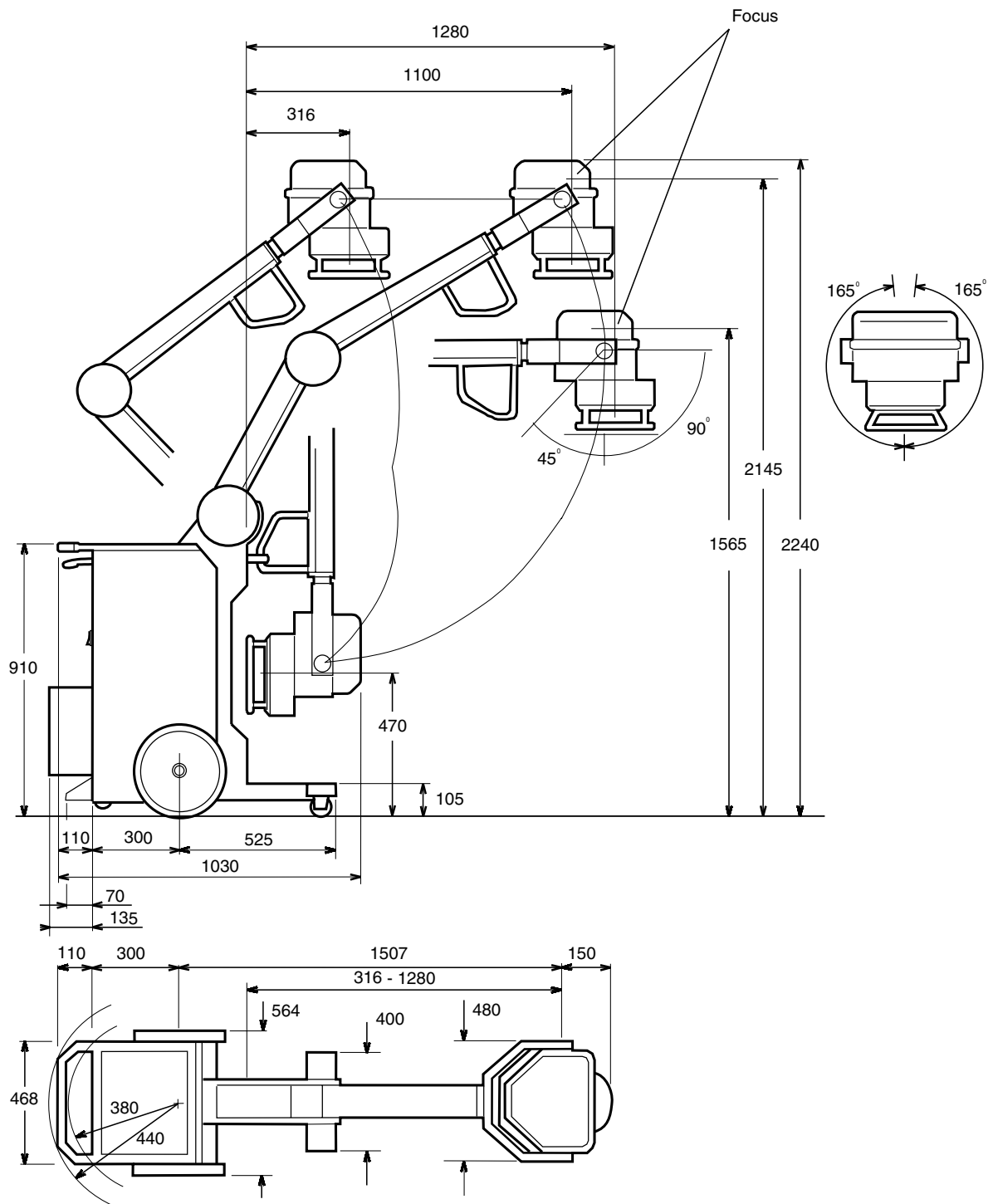

- ◆ Length of the power cable 6 m (with cable winder)
- Length of the power cable 8 m (without cable winder)

This page intentionally left blank.

Electrical data ♦

—	Power supply	R _i	Internal fuse	Power consumption
MOBILETT Plus/ -E/ -M	200 - 240 V ± 10 % 50/60 Hz ± 1 %	max. 2.5 Ohm	10 A at 230 V	max. 1.3 kVA
	100 - 130 V ± 10 % 50/60 Hz ± 1 %	max. 0.6 Ohm	15 A at 110 V	

Weight and heat dissipation

—	Weight [kg]	Heat dissipation [W]
MOBILETT Plus	approx. 210 kg	< 100 W
MOBILETT Plus E	approx. 200 kg	< 100 W
MOBILETT Plus M	approx. 245 kg	< 100 W

Environmental conditions ♦

MOBILETT Plus/ -E/ -M	Operation	Transport / Storage
permissible ambient temperature	+ 10° C ... + 40° C	- 20° C ... + 70° C
permissible relative humidity	30 % ... 75 %	10 % ... 100 % *1

* 1 Maximum 15 weeks.

Packing and transport routes

largest crate	MOBILETT Plus/ -E/ -M	L 1360 x W 780 x H 1860 mm
Weight	MOBILETT Plus	approx. 295 kg with packing
	MOBILETT Plus E	approx. 285 kg with packing
	MOBILETT Plus M	approx. 330 kg with packing

Surface colors ◆

Main color	White textured, Med surface No. 4146 similar to RAL gray - white 9002
Combination color	neutral gray textured, Med surface No. 4426 similar to RAL dust gray 7037
	medical blue, Med surface No. 564 E
Info on spray cans and paint sticks can be found in the projectmanager-handbook (PMH) capital 1	

NOTICE

The MOBILETT Plus/ -E/ -M was not designed for use in the presence of flammable anesthetics (i.e. it is not "anesthetics proof" ["AP"]).

Chapter	Page	Change
0-4		Revision level changed from 03 to 04
1	1-1 and 1-2	General notes updated
2	2-1	Text updated
3	3-1 and 3-2	Technical data updated
4	4-1	Project management, chapter removed
4	4-1	Changes to previous version updated

This page intentionally left blank.